VIVA

ACTIVITY | SPORTS | LIFESTYLE

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MGN1

VERSION: B-BU

www.AsVIVA.de

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WARRANTY TERMS

Important!

Only the original packaging ensures adequate transport protection of the device and should therefore be kept at least in the statutory return period.

For our equipment we provide the legal warranty of 2 years.

For spare parts and repair orders, please contact your dealer or directly to AsVIVA, indicating the customer number or invoice number. Our service is available 24 hours.

Our service address: <https://www.asviva.co.uk/contact>

To ensure that our service staff can help you optimally, please consider the following procedure:

1. Is a specific part damaged

Please tell us the exact part number using the parts list.

2. Is there a fault with the unit?

Should you have the possibility of a video recording e.g. mobile phone camera, pickup the error and send us the video or picture material.

As-STORES GmbH
Kohnacker 9a
D-41542 Dormagen

Homepage: <http://www.asviva.de>


IMPORTANT INFOABOUT YOUR PRODUCT

Thank you for being our valued customer. We are grateful for the confidence you have placed in us and hope we met your expectations. In order to deliver an impeccable product, all our products are subjected to a comprehensive quality control.

If nevertheless an item is defective and needs to be repaired, you do not like the product or you have any suggestions about it, please do not hesitate to contact our support team at service@asviva.de. We will find a solution for your request together.

WHAT NEEDS TO BE CONSIDERED?

- Our company issues a right of revocation for a period of 30 days. Please keep the shipping packaging during this period.
- Please note that during the right of revocation-period the products should only be tested and tried out. Using products over a period of 30 days or cover a mileage over 5km (for example bikes or indoor cycles) does not fall in the usual test phase.
- Your order was delivered packed and free curbside, therefore in case of revocation or defect it will be picked up again packaged and free curbside.
- If you have no more packaging, please contact us and after consultation we can send you one (Notice: reception of the packaging and collection of the device do not occur on the same day.
- The repair of a product takes place within our company. A local service is not possible.

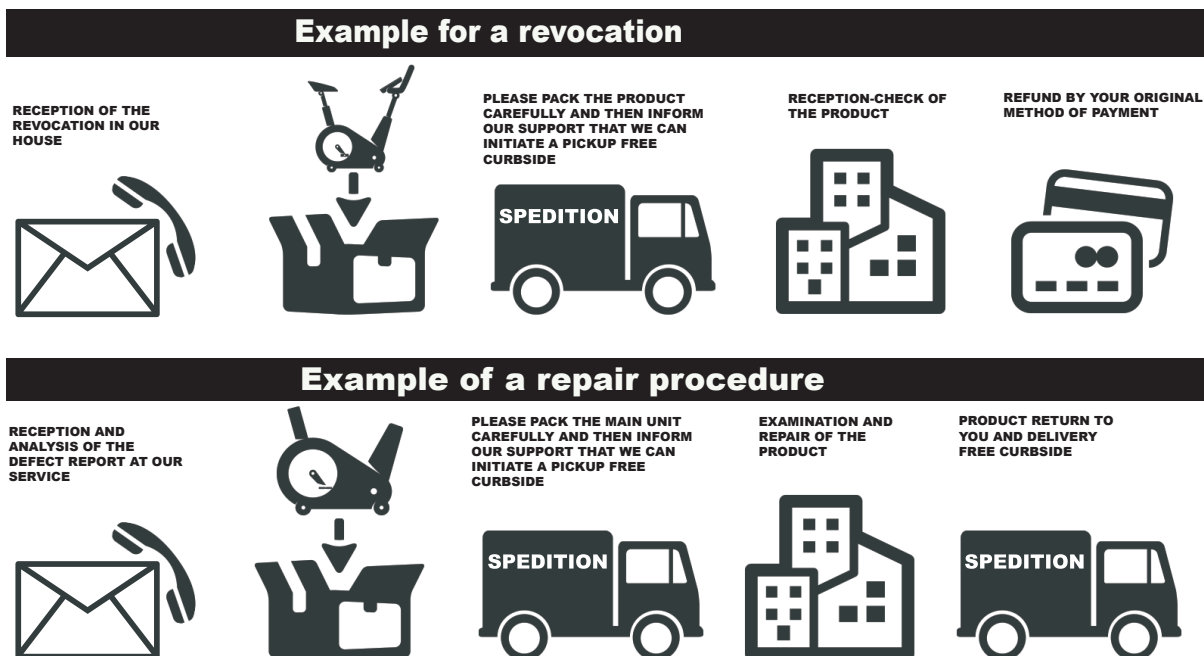IMPORTANT SAFETY INSTRUCTIONS

- The safety level of the device can only be held if it is regularly tested for damage and wear. Particular attention must be paid to the end caps / castors of the feet, the pedals, the handles and the moving elements.
- If you pass this device or use it by another person, make sure that the user is aware of the contents of this manual.
- This unit can only be used by one person to exercise.
- Check that all screws, bolts and other connections are firmly seated before using them for the first time and then at regular intervals.
- Remove all sharp objects from the vicinity of the unit before starting the workout.
- Exercise only on the device if it works properly.
- Defective parts must be exchanged immediately and / or the appliance cannot be used until the maintenance.
- Parents and other supervisors should be aware of their responsibilities because, due to the child's natural playfulness and the ability to experiment, the children will be faced with situations for which the training equipment is not built.
- When you have children trained on the device, you should consider their mental and physical development and, above all, their temperament. You should supervise the children, if necessary, and should point out the correct use of the appliance. As a toy, the training equipment is by no means suitable.
- When installing the unit, please ensure that there is sufficient clearance in each direction.
- Children should not be near the unit unattended, as they can use the unit in a kind and fashion due to their natural play, for which it is not intended and can cause accidents.
- Please note that incorrect and excessive training can endanger your health.
- Please note that levers and other adjustment mechanisms do not extend into the range of motion during the exercise and interfere with the flow.
- When setting up the appliance, make sure that the appliance is stable and possibly ground unevenness are compensated for.
- Always wear training clothes and shoes that are suitable for a fitness training when you are working on the device. The clothing must be designed in such a way that it cannot get stuck anywhere during the workout due to its shape. The shoes should be chosen according to the training device, give the foot a firm hold and have a non-slip sole.
- In general, you should consult your doctor before taking a training course. He can provide you with concrete information about the stress intensity that is appropriate for you and give you tips on training and nutrition.

IMPORTANT SAFETY INSTRUCTIONS

The massage gun is only suitable for adults and should not be used as a toy.

- The massage gun is intended exclusively for private use and is not suitable for commercial use.
- The massage gun is not suitable for therapeutic or medical purposes.
- If you pass this device on or let someone else use it, make sure that he or she is familiar with the contents of these operating instructions.
- Only use it on a dry surface of the body above clothing by gently squeezing it and moving it over the body for about 60 seconds per region.
- Do not use on the head or on hard or bony parts of the body.
- Bruising can occur regardless of the control settings or the pressure applied. Check the treated areas frequently and stop immediately at the first sign of pain or discomfort.
- Keep fingers, hair and other body parts away from the shaft and the back of the massage heads.
- Do not immerse the massage gun in water and do not allow water to enter the ventilation openings.
- Only charge the device with the charger supplied.
- Before each use, check the massage pistol to ensure that it is in perfect condition.
- If the device is visibly damaged, the power cord or plug is damaged, the massage gun must no longer be used.
- Keep the massage gun away from fire, direct sunlight and hot surfaces.

PRODUCT OVERVIEW

- powerful brushless high-torque motor with quiet glide technology
 - rechargeable lithium-Ion battery
- multi-Speed setting for optimal pressure
 - cordless and lightweight (2,4 lbs)
 - interchangeable massage heads

Charge

110 - 240V 50/60 Hz - 5A - 26VDC - 1A

Vibration level

20.0 Hz - 53,3 Hz

PRODUCT OVERVIEW

The massage gun has 7 different heads that can be easily interchanged depending on the area to be massaged.

Fork

For medium-sized muscle groups
(Cervical spine, calves)

Bullet

For the precise treatment of muscles, trigger points and small muscle areas (feet, wrists)

Flat

For general use on large muscle groups

Round

For large muscle groups (thighs, glutes)

Thumb

Suitable for sensitive muscles

Wedge Kopf

Suitable for massaging shoulders, arms, back of thighs and calves

Heat

Suitable for deep relaxation of the whole body muscles

Speed Grade	Speed	Frequency
Small Level	1-129 Level	30-41Hz
Middle Level	130 Level	42Hz
High Level	131-260 Level	43-53Hz


OPERATION

Choose the right massage head for your intended application.

The massage heads are simply pressed onto the message gun until they click into place.
It is not necessary to turn the massage heads.

Switch on the device by turning the main switch on the underside of the handle to the "ON" position.

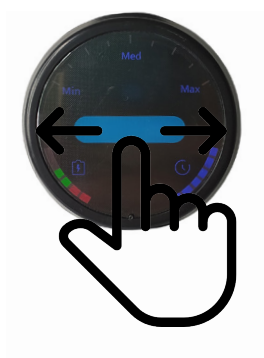By swiping the cruise control
on the right you can increase the intensity levels.

Swiping left decreases the intensity.

- Always exert only light pressure on the massage area with the massage gun.
- Do not massage any point for more than 1 minute.
- Do not use high speeds and / or high pressure on sensitive areas
- Repeated use at too high a speed and / or too much pressure can cause abrasions.
- If you feel pain, stop using it immediately.

After finishing the massage, slide the main switch to the "OFF" position to switch off the massage gun.


BATTERY HANDLING

The green bars on the battery indicator indicate that the battery is sufficiently charged.

The green bars decrease accordingly during operation.

If the red bars light up, the battery must be charged immediately:

- Make sure that the charging socket and the charger are dry.
- Make sure that the massage gun is switched off (main switch to "OFF").
- Connect the charging cable to the charging socket.
- Connect the mains plug of the charger to a suitable socket (220-240 V / 50-60 Hz).
- The LED on the charger glows red while charging.
- If the LED changes to green, the battery is fully charged. You can now finish the loading process.
- Disconnect the charger from the power supply and from the massage gun.

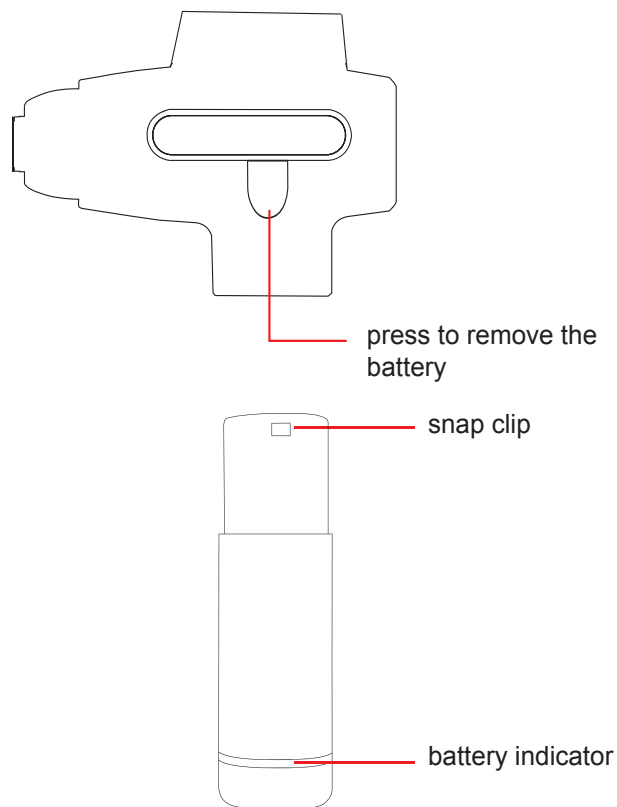The battery can be replaced as required.

To do this, press the button to remove the battery and at the same time pull the battery down and out.

To reinsert, simply plug in the battery until the snap-in clip has clicked into place.

MAINTENANCE AND CLEANING

If you want to care for your device, do not use any aggressive cleaning agents for cleaning, but remove dirt and dust with a soft, possibly slightly damp cloth.

However, avoid contact of liquids of any kind with the device or the interior of the device. This can lead to significant functional impairments.

Information on the battery regulation:

In connection with the sale of batteries or rechargeable batteries or with the delivery of devices that contain batteries or rechargeable batteries, we are obliged to inform you of the following in accordance with the Battery Ordinance: Batteries must not be disposed of with household waste. As the end user, you are legally obliged to return used batteries. You can return batteries free of charge after use at the point of sale or in the immediate vicinity (e.g. in communal collection points or in stores). You can also return batteries to us by mail. Batteries or rechargeable batteries that contain harmful substances are marked with the symbol of a crossed-out garbage can, like the symbol opposite..

The chemical name of the pollutant is located near the garbage can symbol. "CD" stands for cadmium, "Pb" for = lead and "Hg" for mercury.

WARRANTY DEVICE

Congratulations on your decision to purchase a device from AsVIVA.

According to the European guarantee law, amended as of January 1, 2002,

You have a legal warranty period of 2 years.

The guarantee period begins with the handover of the device by the specialist dealer.

In order to prove the date of purchase or handing over, please keep the purchase documents such as the invoice and the cash note for the duration of the warranty period.

In addition to the basic two-year warranty period, there is a statutory

So-called reversal of the burden of proof for a period of 6 months after delivery of the product. This means that the seller is responsible for proving, if an error occurs,

That the defect was not present at the time of handing over and was not caused by improper use, function-related wear or abuse.

1. The warranty does not cover:

- all parts of the device which are subject to wear due to function, in so far as they are not production or material defects.
- Damage caused by improper or faulty maintenance and repairs, alterations or replacement of parts which have not been carried out by a qualified person.
- Damage caused by improper use and force majeure.

2. Authorized warranty claims are provided if:

- The claim or defect at the time of handing over to the customer already existed.
- No natural or function-related wear Cause of the wear or alteration of the goods (see list of wear parts)!
- The damage or defect is not caused by the fact that the device was not used as intended.

If you have questions about your product, our staff will help you
Service department.

We wish you lots of fun and success with your new device!

Your ASVIVA Tteam